


# WAIST MASSAGE

heating / air pressure / vibration /

Red light / Timing massage /

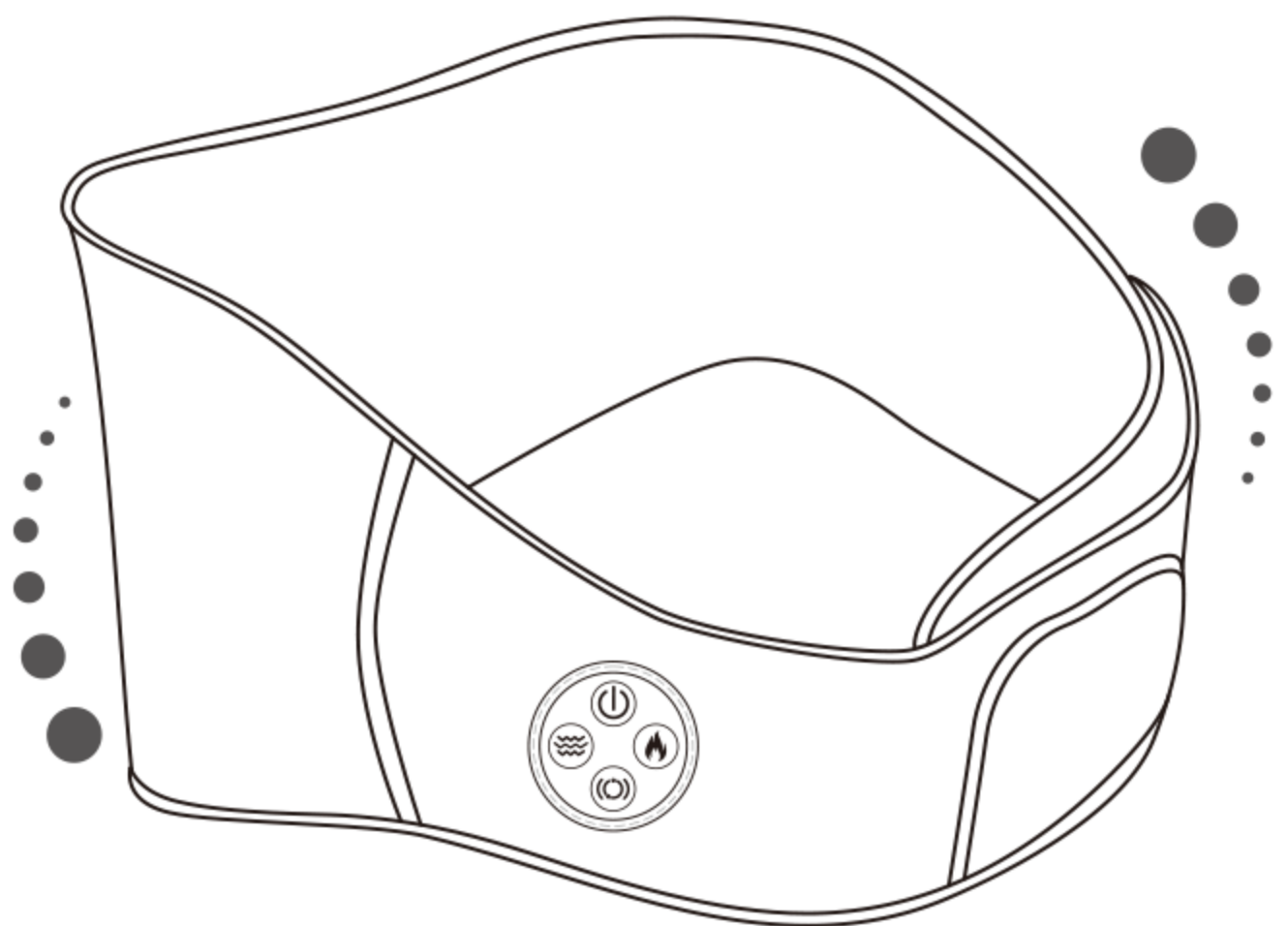

operation instruction and user manual

## PRODUCT INSTRUCTION AND FUNCTION:


- ①USB port connect 5V 2A batttery bank or adaptor.
- ②easily carry and enjoy massage.
- ③product surface is good looking and fashion, leading trend.
- ④multi-functional massage. Red wavelength 650nm effectively alleviate discomfort
- ⑤Maxium rated power 8W and working long time.




## USE PRODUCT METHOD

### 1、pre-use preparation

(1)take it out

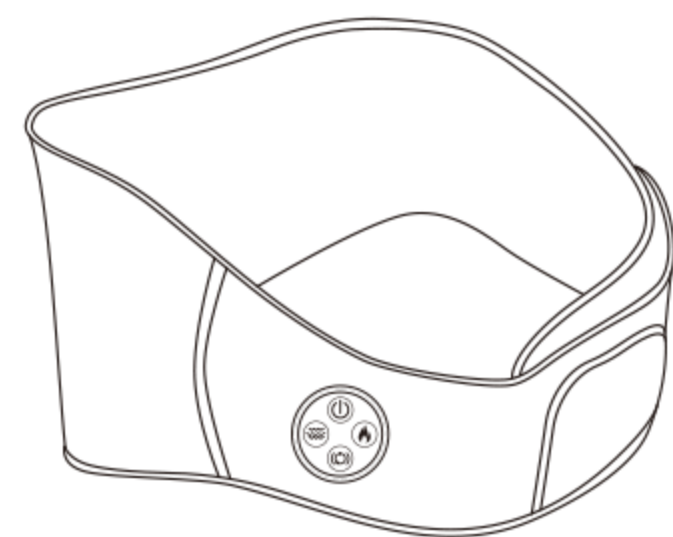


(2)connect battery bank or adaptor,press the power key 2 sec

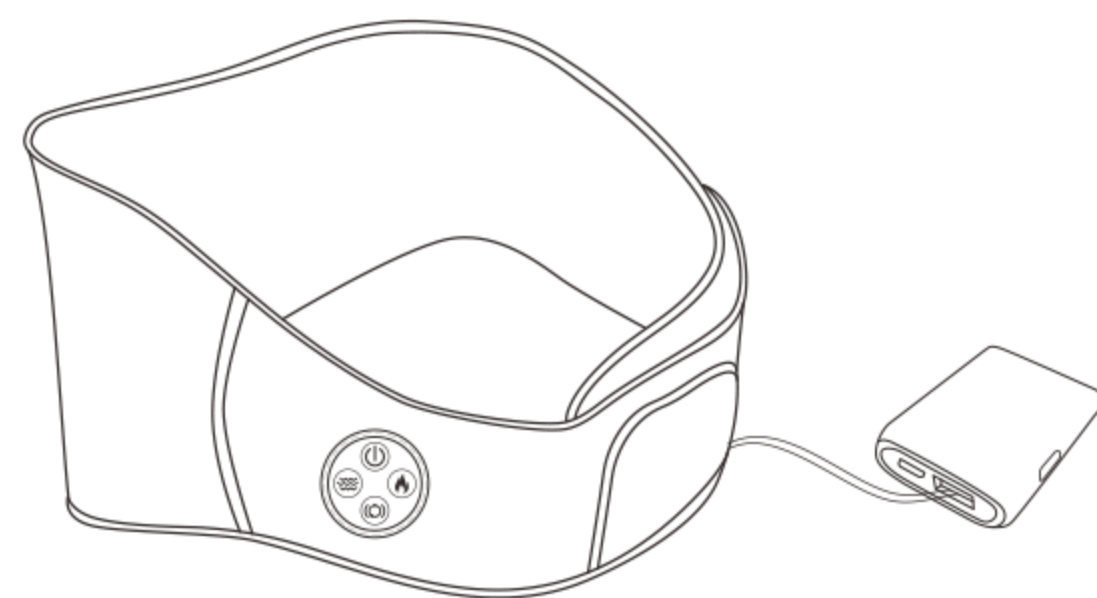


## 2、operation instruction (below picture)

(1) open magic sticker,bundle it to comfortable position

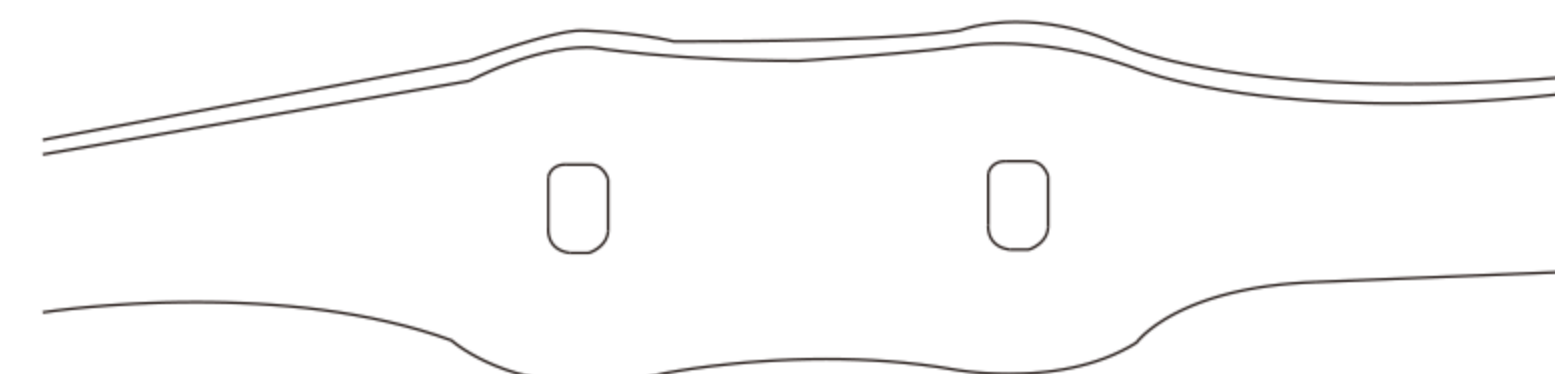

(2) use magic sticker to adjust right position

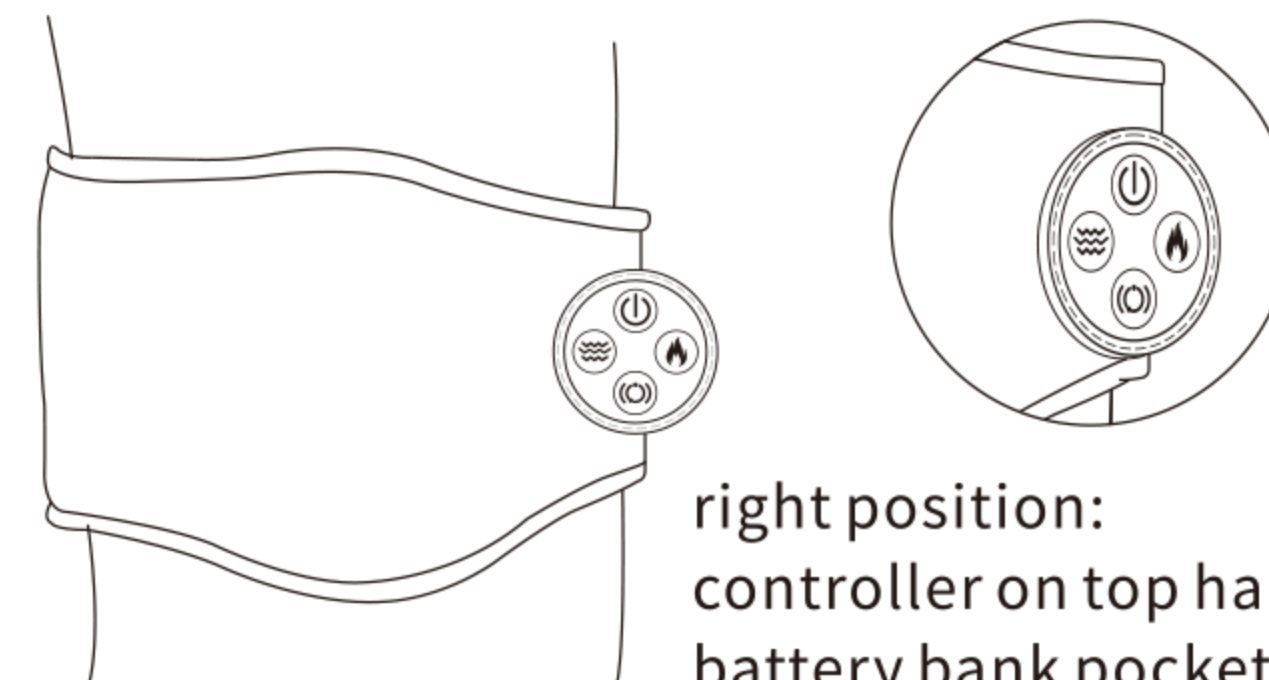

right position:  
controller on top half,  
battery bank pocket on lower half.

(3)press power key 2 sec.all function working,green light open


(4)press heating key:low/middle/high/off(light color:green/blue/red)

(5)press air pressure:Low/High/off(light. color:green)

(6)press vibration:weak/pulse/steady/off(light color:green/blue/red)

(7)automatic off when continuous working 15 minutes

## 3、controller instruction

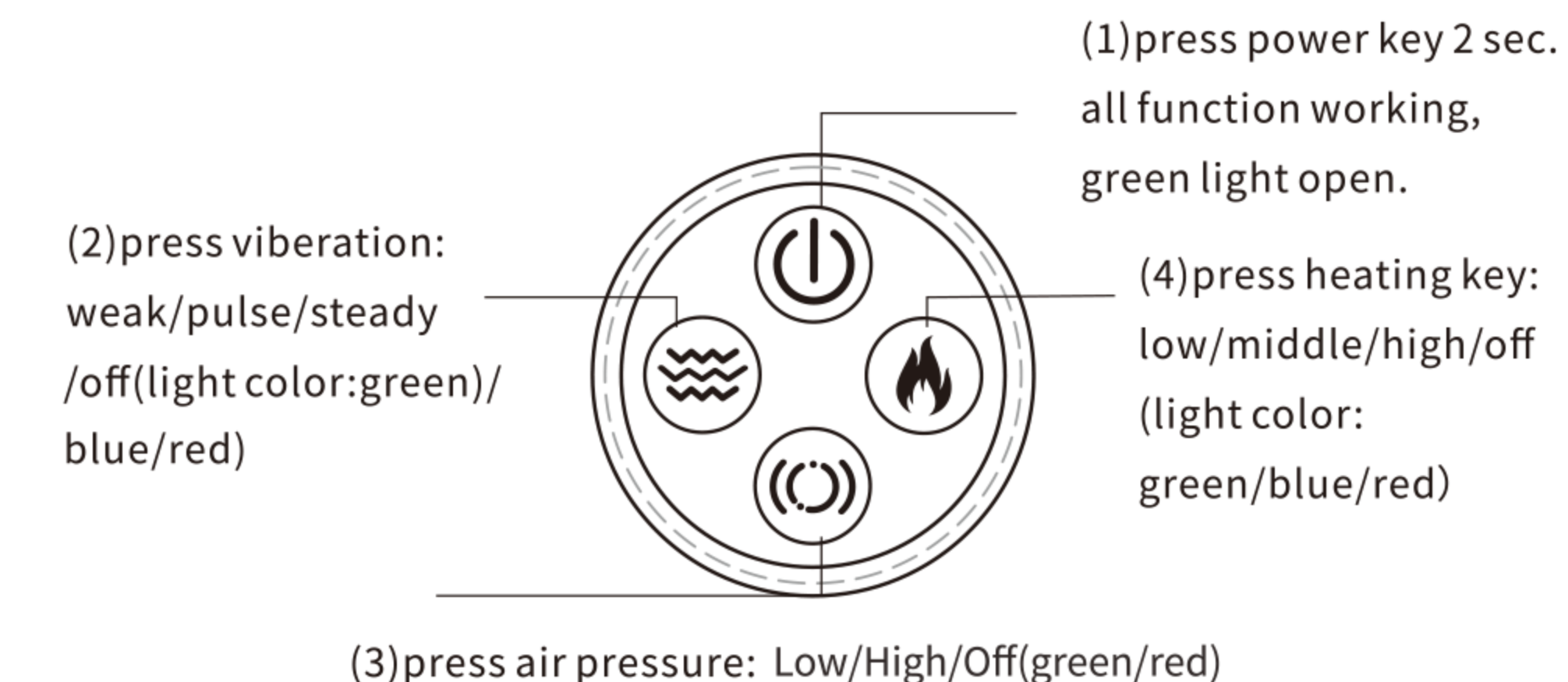

## 4、precautions!

(1)Be careful to use the highest heating for long periods of time to avoid burning skin.

(2)prohibited use product when sleeping

(3)can wipe surface of product,prohibited into water to clean

## 5、disable crowds:

(1) skin with thermal allergies

(2)unable take care of themselves

(3) doctor warning